





ElevenLabs Integration for Mind Express

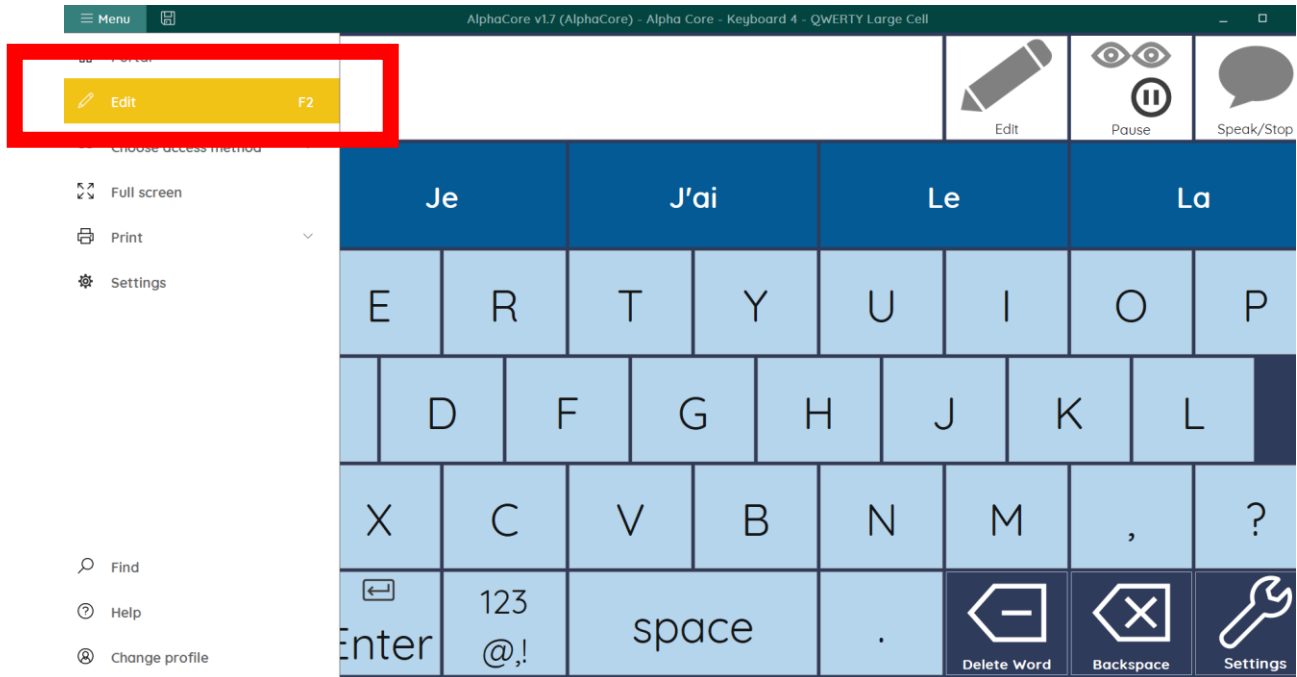


1. Open the main keyboard in Mind Express, navigate to the menu screen and click Edit

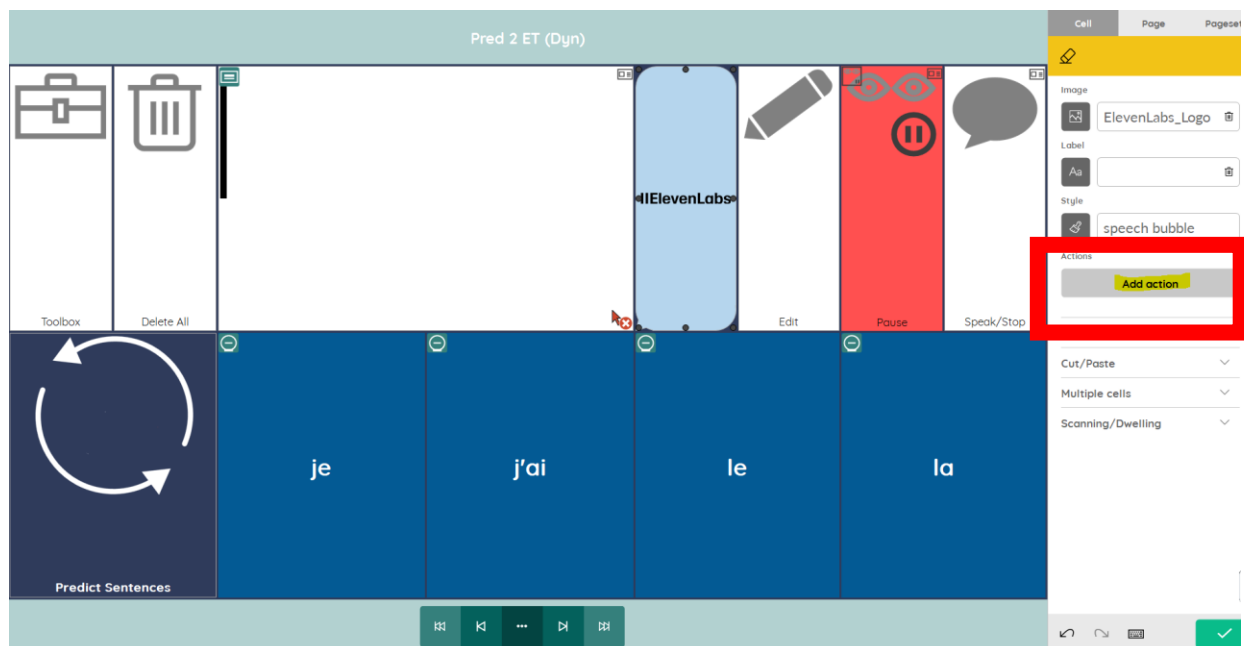


2. Add a label for you button and then click Add Action

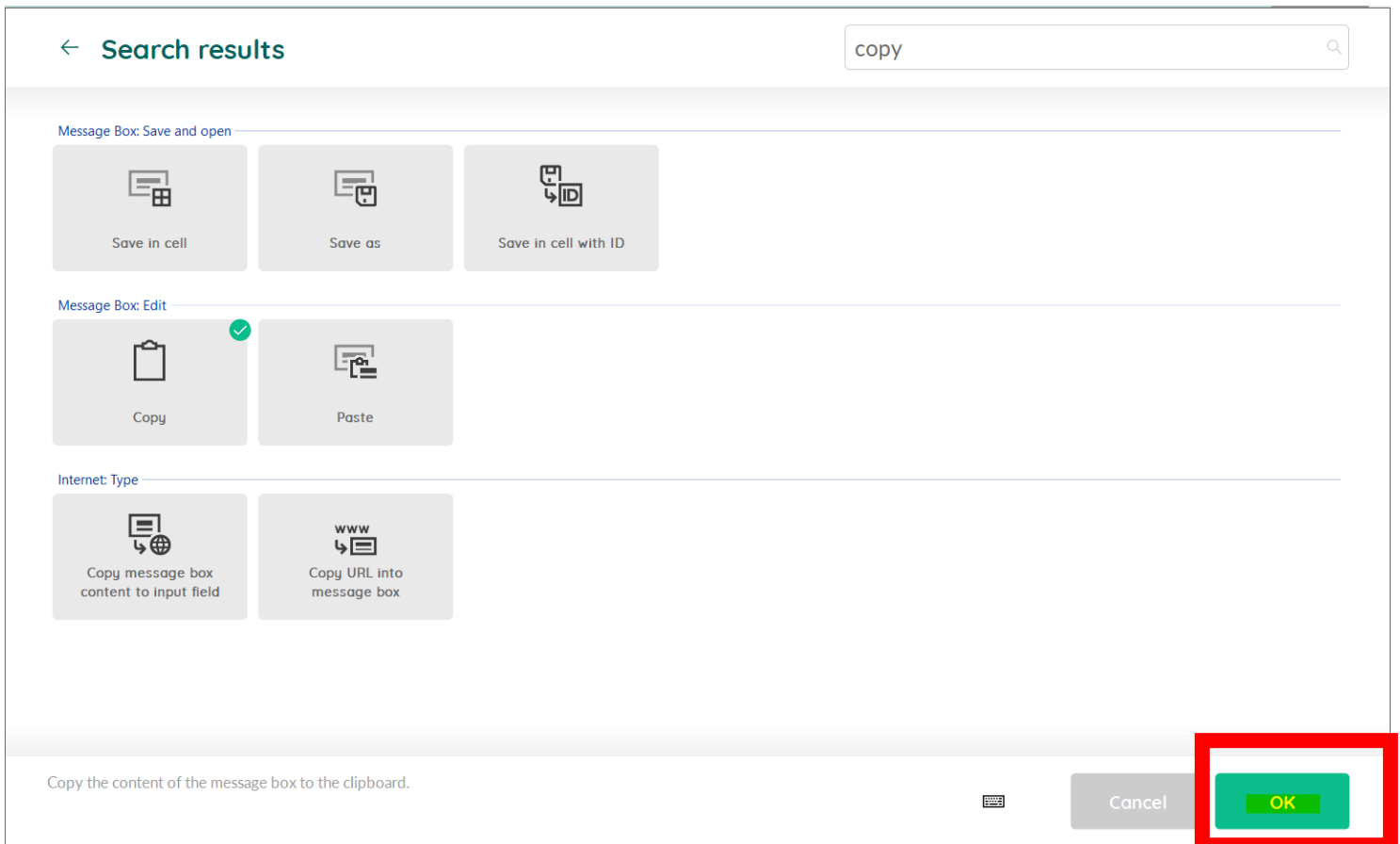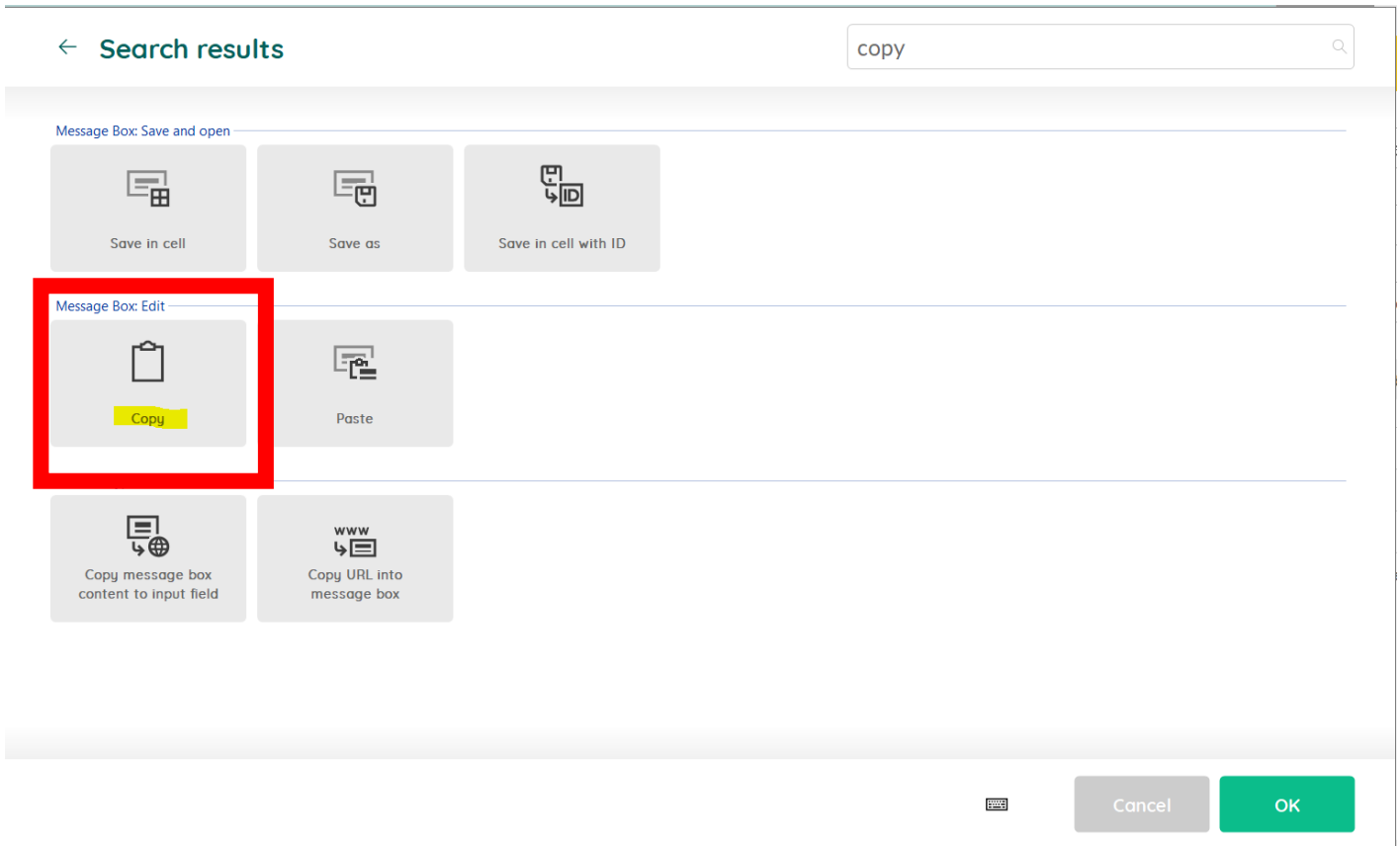


3. Search for Copy and then click OK

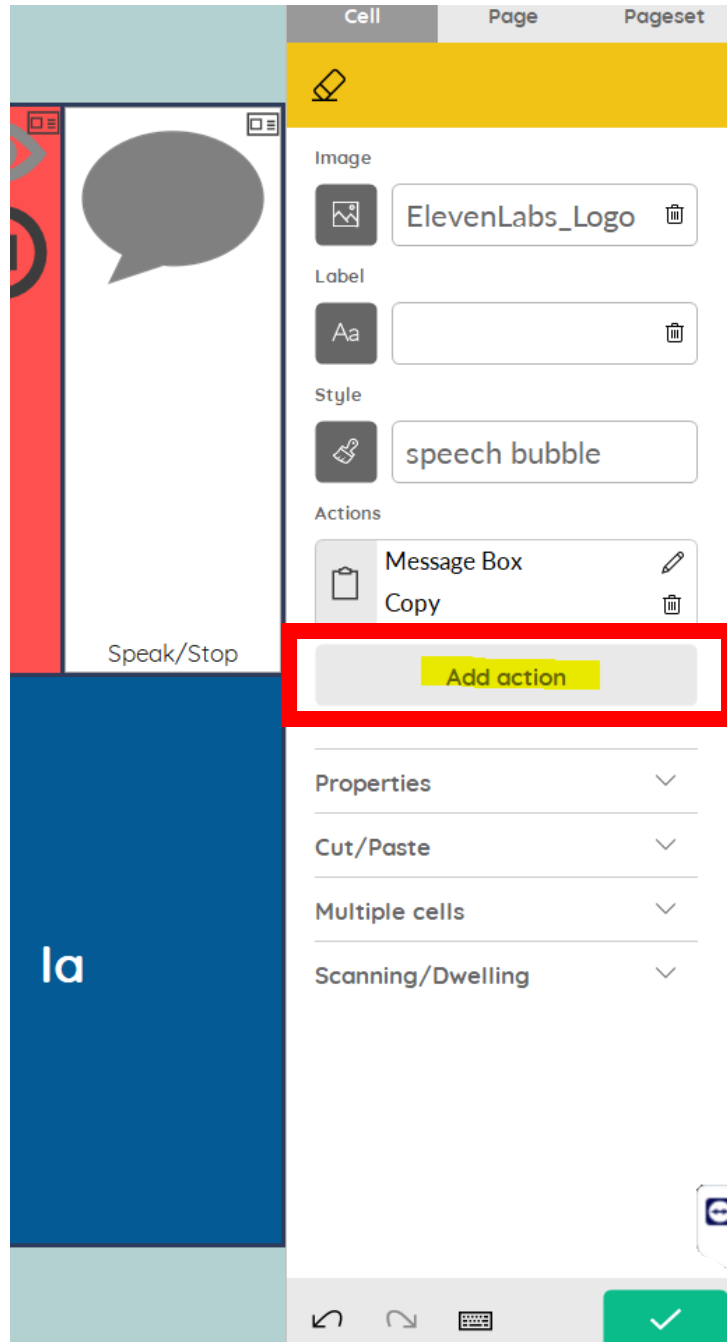


4. Click on Add Action

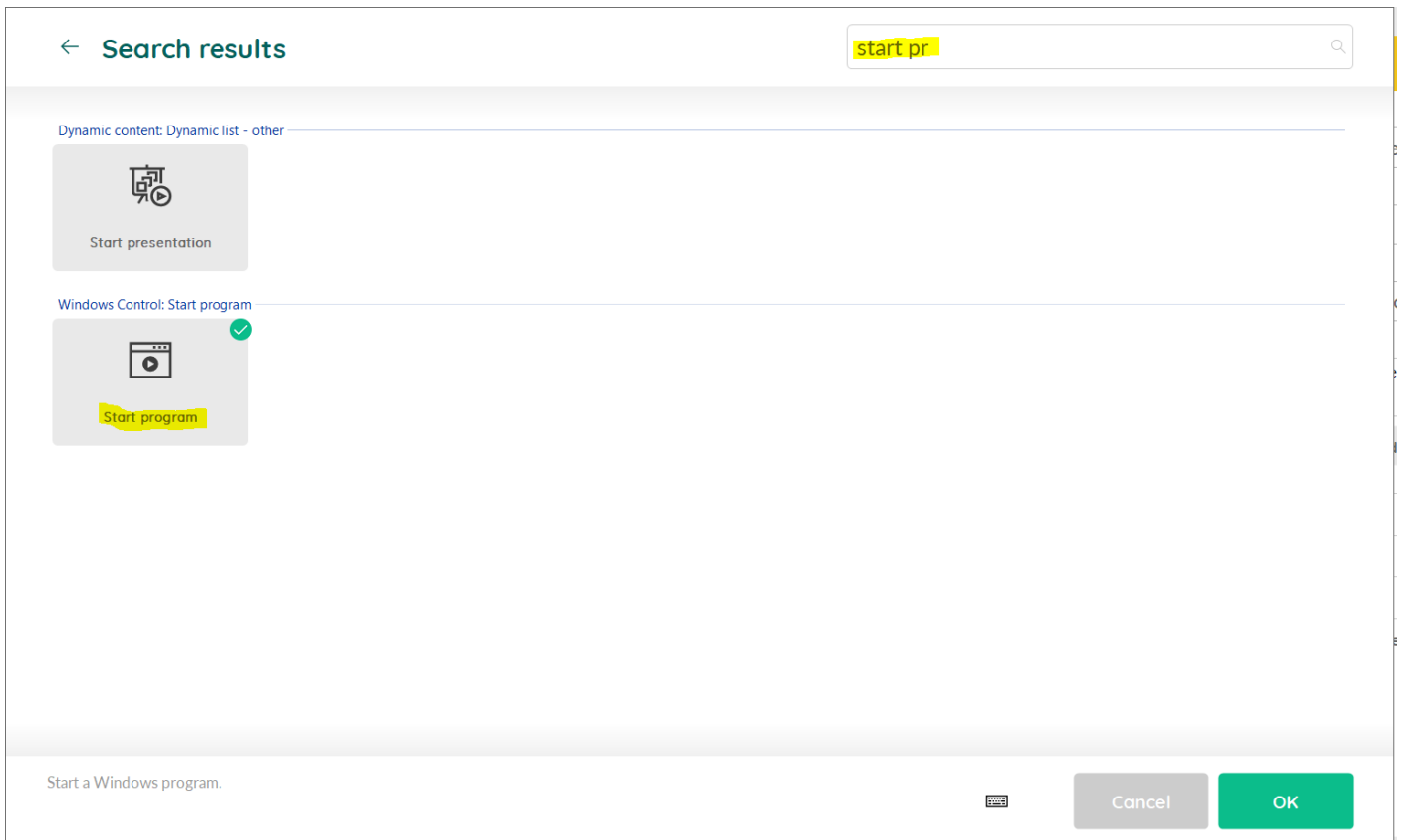

5. Search for Start Program, select and then click OK

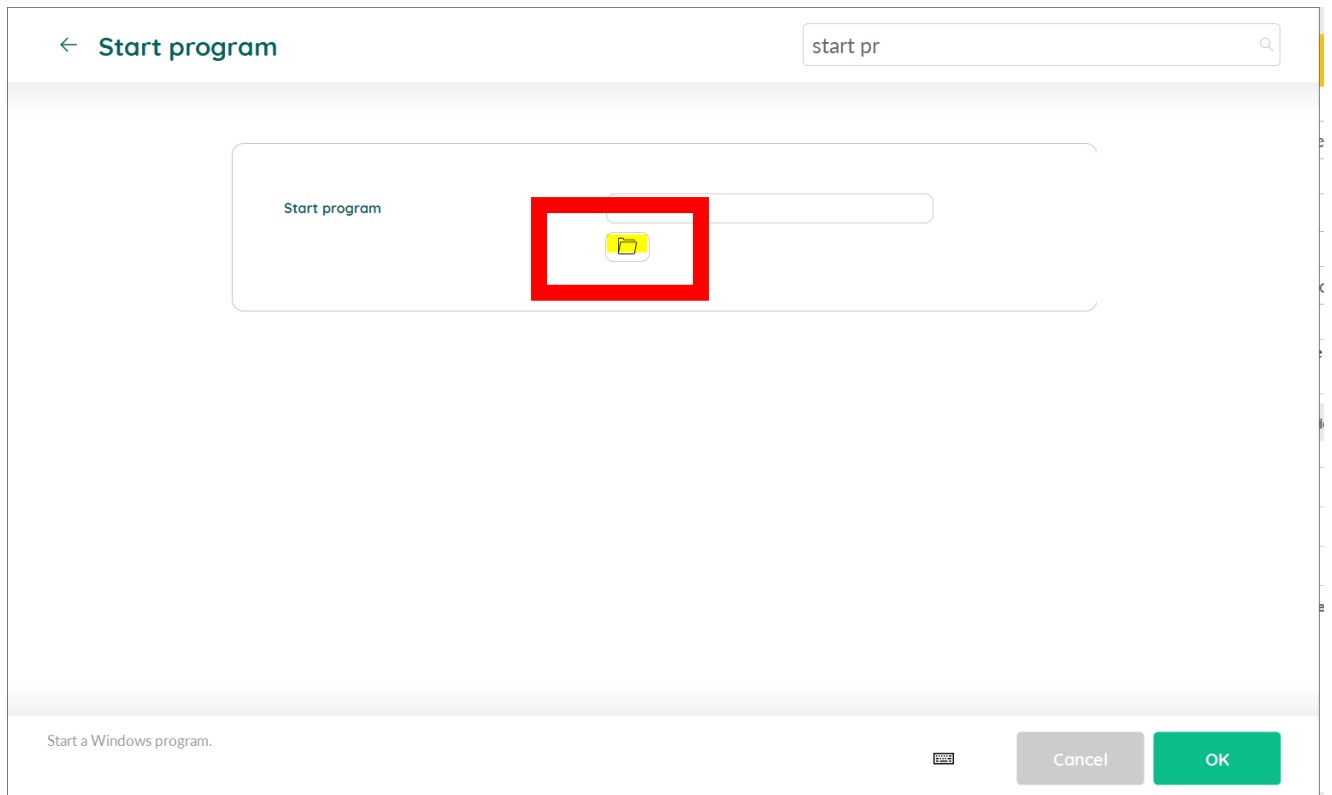

6. Click on the folder icon

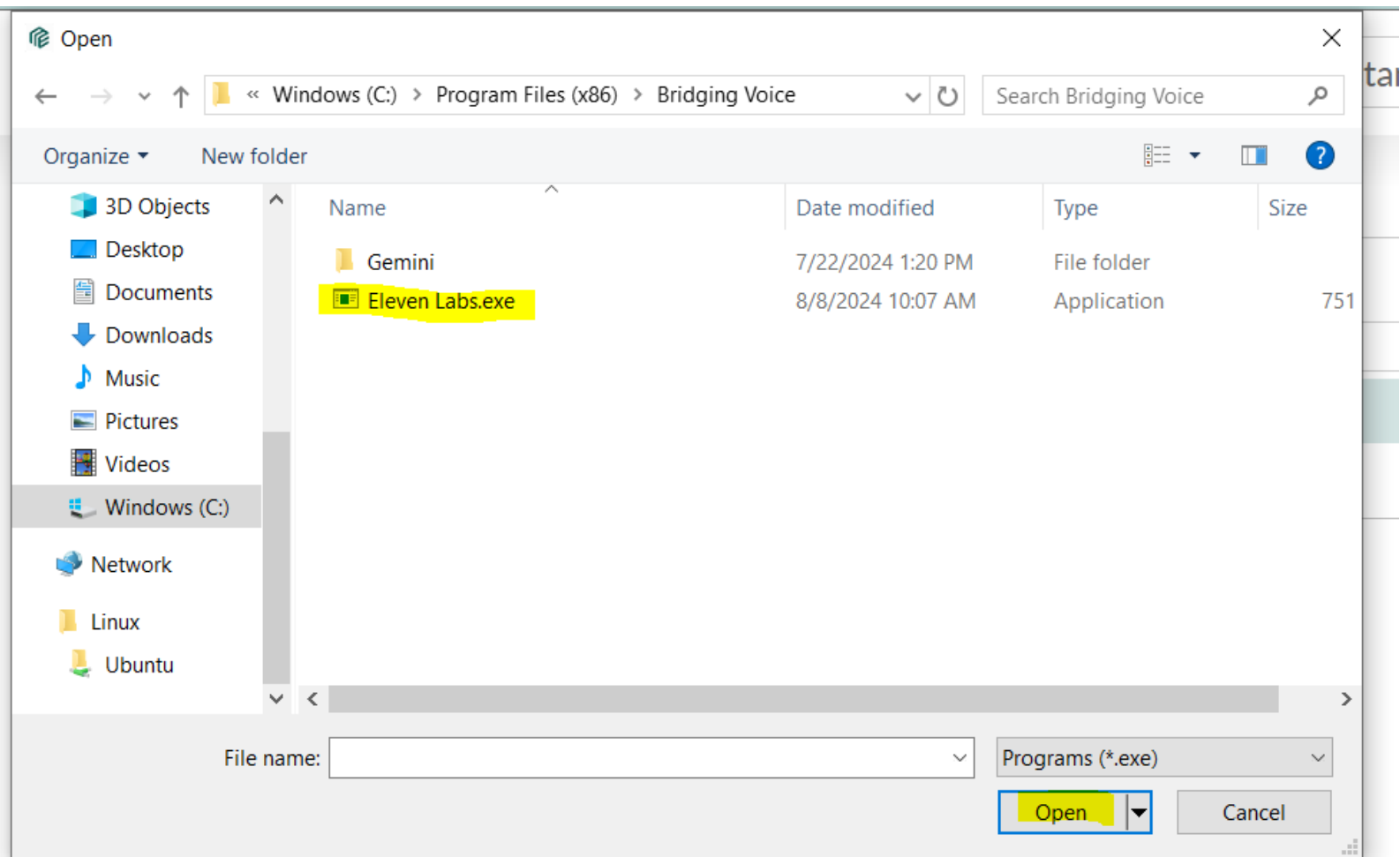


7. Navigate to “C:\Program Files(x86)\BridgingVoice and select ElevenLabs.exe. Then click Open

Start program


oice\Eleven Labs.exe" Xb7hH8MSUJpSbSDYk0k2




8. Using the right arrow, move your cursor to the end of the Start program box and add a space
9. Open the Word or Notepad doc where you stored the ElevenLabs API KEY and Voice ID
10. Highlight and Copy (Ctrl+c) the Voice ID first
11. Back in Mind Express Paste the Voice ID into the "Start Program" (leaving one space after the file name)
12. Go back to the Notepad or Word doc and copy the API Key.
13. Back in Mind Express, paste the API Key into the "Start Program" after the Voice ID with a single space between.

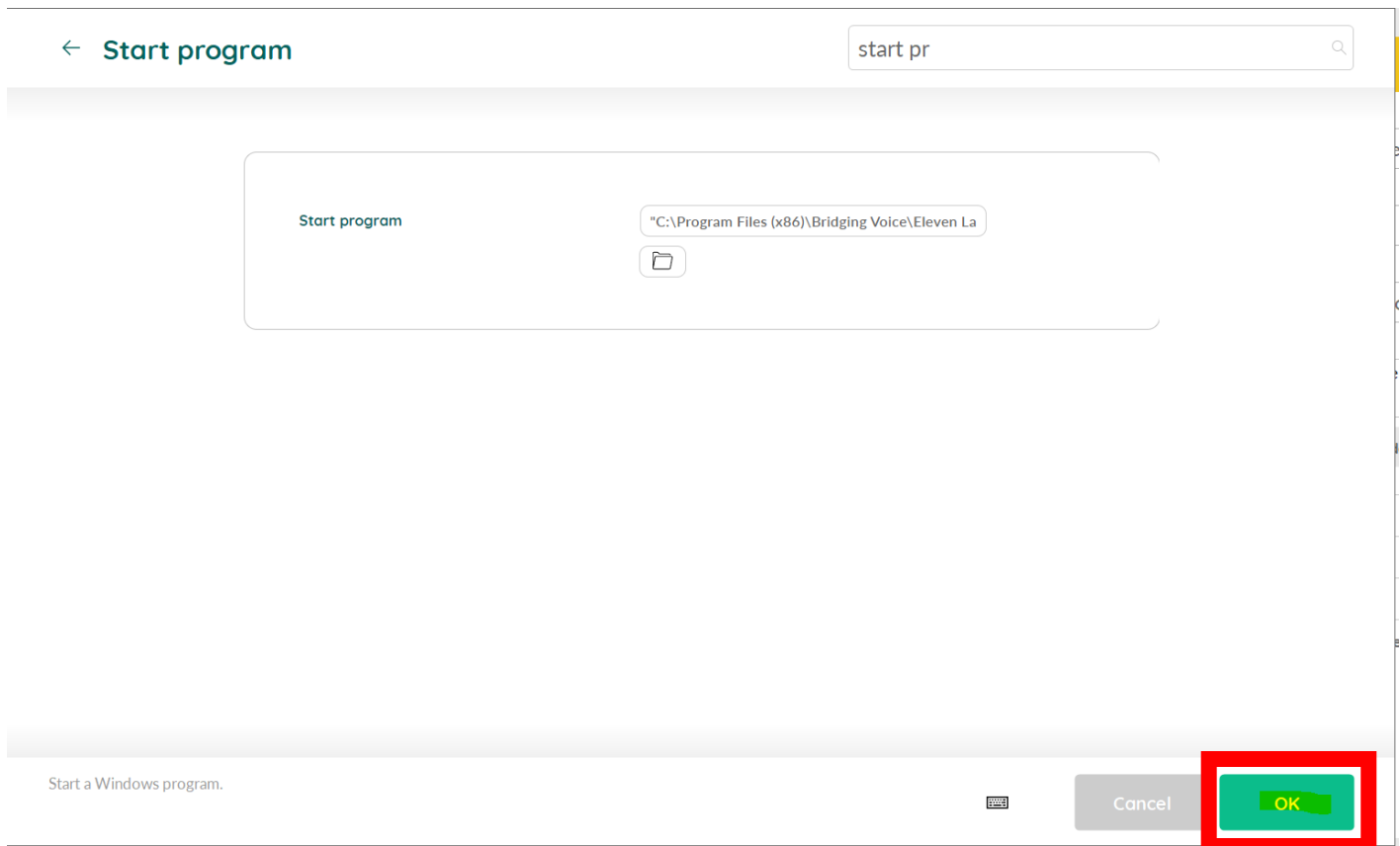

14. Press OK

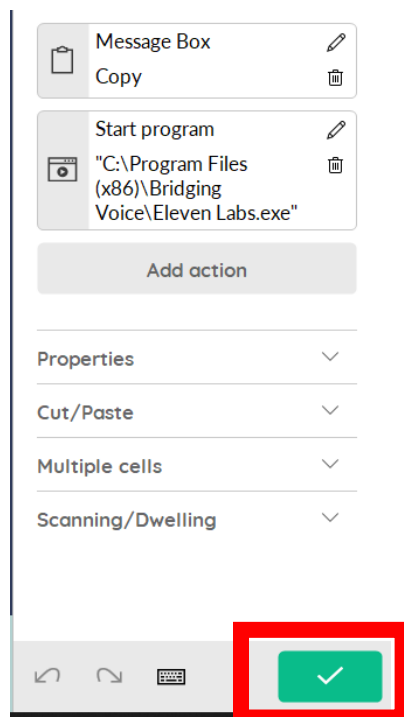

15. Click the Green Check Box in the bottom right to save the button.

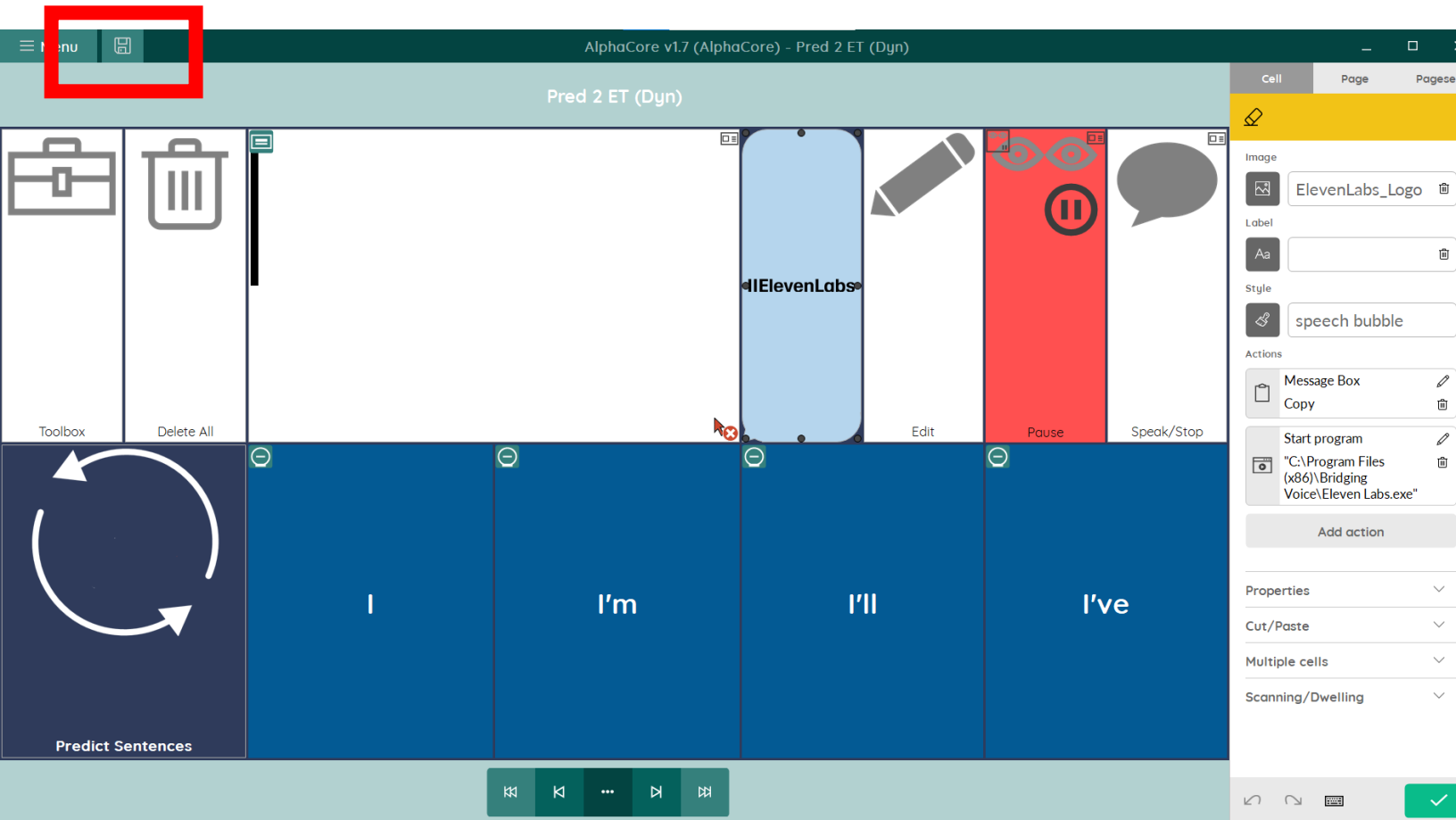

16. Click the little Save icon (floppy disk) next to Menu in the upper left corner.